

Tranqwil



TT850 Eletronic Toothbrush User Manual

TPC-02S


T

For all comments, questions or complaints,
please reach out to us.

Email: hello@tranqwilhome.com
Website: Tranqwilhome.com

Product Introduction

Components

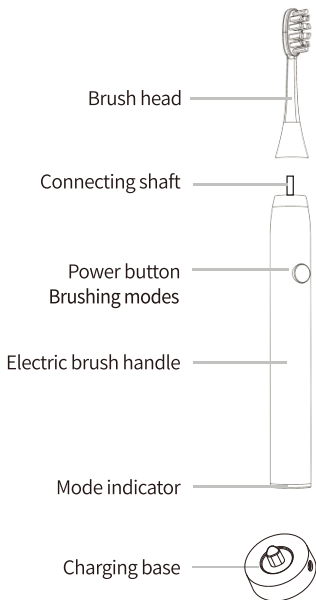

Power Button/Brushing Modes


Turn on: Press the power button to turn on the toothbrush.

Mode switch: After turn on the toothbrush, press the power button again to switch between the 5 brushing mode; after the 2-minute cycle completes, the toothbrush will turn off.

Turn off : Press the power button 5 seconds after a mode being selected, the toothbrush will turn off. The toothbrush will memory the mode you used last time.

5 Modes : The electric toothbrush has five modes and its indicator lights are the followings:

- Low Speed Cleaning Mode : Light blue light
- Medium speed Cleaning Mode: Blue light
- High speed Cleanig Mode: Red light
- Whitening Mode: White light
- Massage Mode: Green Light

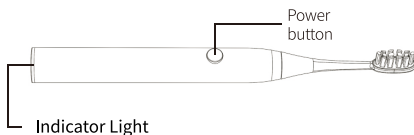

How to Brush

Tranquil TT850 electric toothbrush completely adapts to the brush method recommended by ADA and cleans cavity effectively. Please follow the advice to achieve the best cleaning result.



Brush head keeps the angle of 45 degrees from teeth to clean out dental plaque gingival sulcus better.




Intelligent alarm with short pause for changing area every 30s



Increasing strength prevents from splash of water, comfortable during oral cleaning


Cleaning and Maintenance

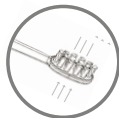

Clean brush head and keep ventilation



Completely water proof, feel free to wash with running water

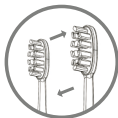

Suggest to change brush head every 3 months

Charging

When the colored indicator lights at bottom glitter during usage please charge it as soon as possible.

Please note:

1. Do not immerse charging base to water or any other liquid.
2. Do not place charging base in high temperature environment for charging, usage or custody.
3. Do not leave it on charging base for a long time for the best use of battery service life.
4. Do not use the product when charging.



- Green indication, 20%-100% power, please feel easy to use.
- Red indication, 0%-10%, please charge as soon as possible.

Using Instructions

Notes:


Do not charge in bathroom, nor place the charger above or near water source, e.g. bathtub, water sink.

Do not replace any product parts by yourself, please contact Tranqwil for part replacement; unauthorized disassembly of product is strictly prohibited.

Cleaning and maintenance

As shown in the following picture, match the direction of toothbrush bristle with that of the button on brush holder, then assembly them.

Since the bottom of brush holder is fixed on travel-size storage box, please take up the holder with your hands holding at the power button position, and pull it out from the groove.

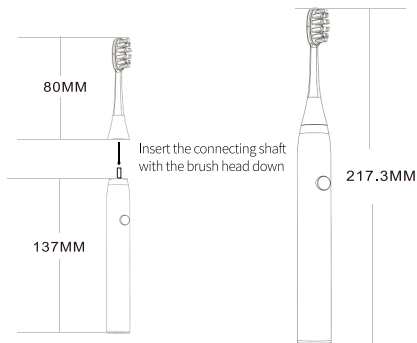

Specifications

Electric brush holder	Product mode no.	Apollo 7
	Rated voltage	3.7V
	Working power	1.3W \pm 0.7W
	Water-proof grade	IPX7
	Charging time	5 hours
Charging base	Input parameter	5V
Date of production	As shown on the label at the back of package	
Place of production	Bao'an District, Shenzhen City, Guangdong	

Troubleshooting

Troubles	Causes	Solutions
Unable to connect brush head with holder correctly	Brush head jack is blocked by foreign matters; brush holder adheres to foreign matters	Clean out foreign matters
Unable to turn on	Low battery or run out of battery	Charge now, and turn on normally
Battery life	The longest time the battery can sustain under standby status when it is fully charged	57 days when it is fully charged
Toothbrush care	Bacteria may grow in old toothbrush head, which may do harm to people when they brush teeth	Clean out foreign matters, replace toothbrush head every 1-2 months
How to clean toothbrush	After brushing teeth, bubbles may be splashed to the holder	Clean by hot water or hot towel
Charging time	Charging under standby status	Full charging needs 5 hours